



- Available with either screw cap for low speeds or leakproof sealing cap for high speeds
- Sealing Cap designed with built-in grips on caps
- Available in Polycarbonate or Polypropylene

Warnings and Guidelines


- May be autoclaved to 121°C at 15 psi for 20 minutes
- Bottles must be filled to 80% capacity for centrifuge work
- Prolonged autoclaving and reuse may cause cracks or discoloration

Roush Life Sciences' **Autofil™** disposable centrifuge bottles have been designed specifically for Life Science researchers to process samples for cell suspensions, cell lysates, tissue homogenates and organ homogenates.

Our centrifuge bottles feature thick walls that hold up to reuse, are available with a leakproof sealing cap to avoid messy rotor cleanup and wide mouths that make filling and pouring liquids easy. Each **Autofil™** centrifuge bottle also features large, clear-molded graduation markings and a frosted area for convenient labeling.

Autofil™ centrifuge bottles, available in either 250 ml or 500 ml, fit a wide range of rotors. Our centrifuge bottles are available in polycarbonate or polypropylene with your choice of a traditional screw cap or a leakproof cap with a sealing silicone O-ring.

- 1 Screw on or leakproof sealing cap with built-in grips (leakproof sealing cap shown)
- 2 Available in autoclavable polycarbonate or chemical-resistant polypropylene
- 3 Crystal clear polycarbonate or opaque polypropylene
- 4 High visibility graduations

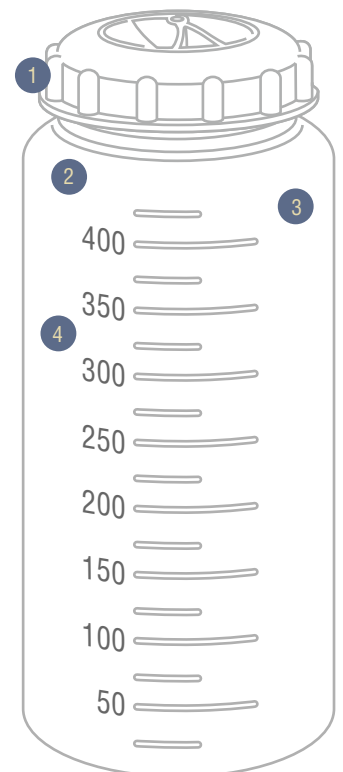



	Item Number	Description	Dimensions	Closure Type	Max. RCF (g) w/cap	Units/Pack	Packs/Case
Polypropylene	1801-RLS	250ml PP Bottle with Screw Cap, nonsterile	2.32 x 5.01 in 59 x 129 mm	Screw	8,250	4	9pk
	1803-RLS	250ml PP Bottle with Seal Cap, nonsterile	2.32 x 5.02 in 59 x 129 mm	Seal	27,500	4	9pk
	1805-RLS	500ml PP Bottle with Screw Cap, nonsterile	2.68 x 6.61 in 59 x 129 mm	Screw	8,250	4	6pk
	1807-RLS	500ml PP Bottle with Seal Cap, nonsterile	2.68 x 6.61 in 59 x 129 mm	Seal	13,500	4	6pk
Polycarbonate	1802-RLS	250ml PC Bottle with Screw Cap, nonsterile	2.32 x 5.01 in 59 x 129 mm	Screw	8,250	4	9pk
	1804-RLS	250ml PC Bottle with Seal Cap, nonsterile	2.32 x 5.02 in 59 x 129 mm	Seal	31,700	4	9pk
	1806-RLS	500ml PC Bottle with Screw Cap, nonsterile	2.68 x 6.61 in 59 x 129 mm	Screw	8,250	4	6pk
	1808-RLS	500ml PC Bottle with Seal Cap, nonsterile	2.68 x 6.61 in 59 x 129 mm	Seal	13,850	4	6pk



Chemical compatibility information is intended as a guide. We recommend that you test bottles for your specific application before use.

Roush Life Sciences, 8E Industrial Way #8, Salem, NH 03079 (603) 540-5778
www.RoushLifeSciences.com