



- **Autofil™** is available as a complete system or as a bottle-top filter
- **Autofil™** design including an ergonomic bottle shape, molded finger grips and patented bottle cap
- Conveniently sold presterilized


The **Autofil™** laboratory filtration system is a disposable, sterile, vacuum operated system used for the preparation of buffers, tissue culture media, microbiological media and other biological fluids.

Autofil™ is available either as a complete system or as a bottle top filter. Both options have an upper filter funnel with a dust cap that contains a 0.22µm (sterilization grade) or 0.45µm (clarification grade) asymmetric PES membrane filter. Our membrane is low protein binding and surfactant free. The complete system includes a Roush Life Sciences **SECUREgrasp™** bottle attached to the filter.



ISO 9001
ISO 13485

- 1 Quick-Connect™ pedestal base or cradle ring allows for hands free filtration
- 2 Forward-facing, large, raised graduation marks for easy reading
- 3 **SECUREgrasp™** bottle cap technology
- 4 Available in 0.22µm or 0.45µm asymmetric membrane filter

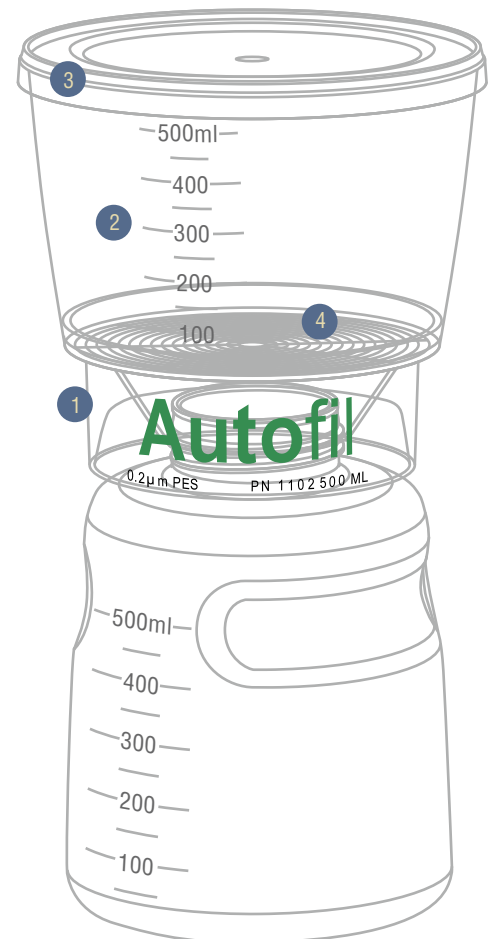



	Item Number	Volume	Pore Size & Membrane	Filter Diameter	Description
Full Assembly	1101-RLS	250	0.22µm PES	66mm	250ml complete system, 0.22µm PES, sterile, 12pk
	1102-RLS	500	0.22µm PES	90mm	500ml complete system, 0.22µm PES, sterile, 12pk
	1103-RLS	1000	0.22µm PES	90mm	1000ml complete system, 0.22µm PES, sterile, 12pk
	1141-RLS	250	0.45µm PES	66mm	250ml complete system, 0.45µm PES, sterile, 12pk
	1142-RLS	500	0.45µm PES	90mm	500ml complete system, 0.45µm PES, sterile, 12pk
	1143-RLS	1000	0.45µm PES	90mm	1000ml complete system, 0.45µm PES, sterile, 12pk
Funnel Only	1151-RLS	250	0.22µm PES	66mm	250ml bottle-top filter, 0.22µm PES, sterile, 12pk
	1152-RLS	500	0.22µm PES	90mm	500ml bottle-top filter, 0.22µm PES, sterile, 12pk
	1153-RLS	1000	0.22µm PES	90mm	1000ml bottle-top filter, 0.22µm PES, sterile, 12pk
	1161-RLS	250	0.45µm PES	66mm	250ml bottle-top filter, 0.45µm PES, sterile, 12pk
	1162-RLS	500	0.45µm PES	90mm	500ml bottle-top filter, 0.45µm PES, sterile, 12pk
	1163-RLS	1000	0.45µm PES	90mm	1000ml bottle-top filter, 0.45µm PES, sterile, 12pk
Hardware	1120-RLS	N/A	N/A	N/A	Cradle ring, 1/pack
	1111-RLS	N/A	N/A	N/A	Pedestal Stand, 1/pack



Chemical compatibility information is intended as a guide. We recommend that you test bottles for your specific application before use.

*US Patent Number: D587580; Indian: 215829; Japanese: 2008-3034
 Roush Life Sciences, 8E Industrial Way #8, Salem, NH 03079 (603) 540-5778
www.RoushLifeSciences.com