

- **VersaCap[™]** technology provides the flexibility of interchangeable top-connectors
- **VersaCap[™]** is offered in industry standard size and also wide-mouth version to fit **EZfil[™]** carboy
- Quick release latches for easy port change
- Leak-proof cap with replaceable gasket
- **VersaCap[™]** is stackable
- Available in a variety of configurations
- **VersaCap[™]**, **VersaBarb[™]** and Spigot are autoclavable
- The RLS Spigot also offers a convenient solution for bottom-connections and can fit both 5/8" and 3/4" hose sizes
- RLS spigot controls two flow rates

Roush Life Sciences has expanded the **EZgrip[™]** carboy features and benefits by offering the **VersaCap[™]**, **VersaBarb[™]** and state-of-the-art Spigot Technology. These three components complement a wide array of applications supported by the **EZgrip[™]** carboy.

VersaCap[™] technology (patent pending) provides the ultimate in top-connection flexibility. The leak-proof cap insert with replaceable gasket enable users to quickly change the top connector of each carboy to match the application needs.

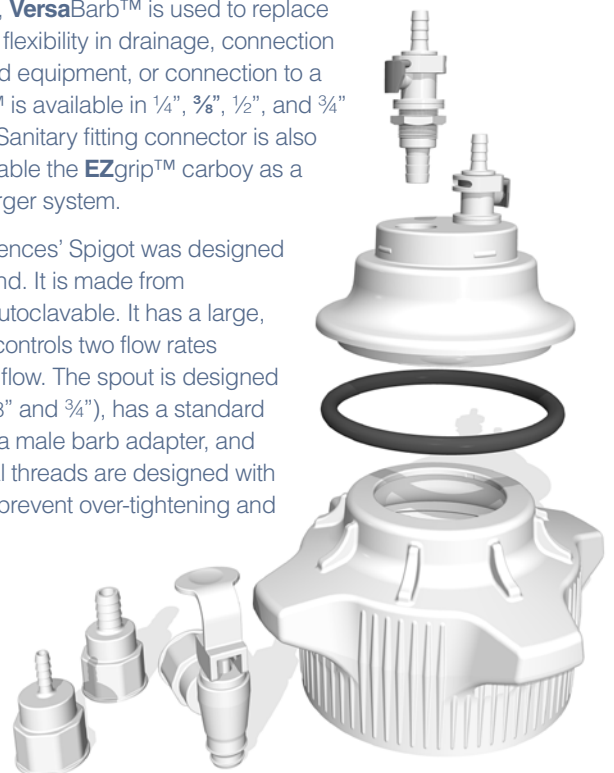
VersaCap[™] technology has large handles with finger grips that allow for easy tightening and opening.

VersaCap[™] is stackable and is offered in both 83mm (compatible with most industry carboys) and a 120mm wide-mouth version that fits exclusively on Roush Life Sciences' 20L and 40L **EZfil[™]** carboy. And to visualize homogenization of solutions, closed adapter inserts with a clear window are available.

Quick-connect options are also a part of the many capabilities of this **VersaCap[™]** technology. This includes quick-connect port combinations of two, three, and four. All quick-connect options are available in 1/4" and 3/8" hose barb sizes, the industry standard and molded hose barb options are available with 1/4" and 1/2" molded hose barbs.


For bottom tubulation, **VersaBarb[™]** is used to replace the spigot and support flexibility in drainage, connection to a pump or automated equipment, or connection to a filling bell. **VersaBarb[™]** is available in 1/4", 3/8", 1/2", and 3/4" hose barb sizes. A 3/4" Sanitary fitting connector is also available in order to enable the **EZgrip[™]** carboy as a supply reservoir to a larger system.

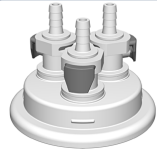
Lastly, Roush Life Sciences' Spigot was designed with the end-user in mind. It is made from polypropylene and is autoclavable. It has a large, comfortable lever that controls two flow rates including a continuous flow. The spout is designed to fit two hose sizes (5/8" and 3/4"), has a standard internal luer taper to fit a male barb adapter, and the spigot body internal threads are designed with a positive hard stop to prevent over-tightening and damage to the seal.







VersaCap™ 80mm


	205-1003-RLS			
	Cap, Closed			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	n/a	n/a	2

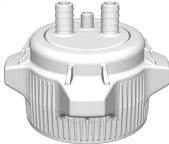
	205-3004-RLS			
	Adapter, 1/4" HB x 3			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	1/4"	3	1

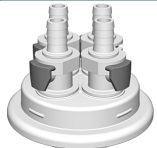
	205-1004-RLS			
	Cap, Open with closed adapter			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	n/a	n/a	2


	205-3005-RLS			
	Adapter, 3/8" HB x 3			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	3/8"	3	1


	205-1007-RLS			
	Cap, Open top with adapter, 1/4" HB x 2 and Vent			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	1/4"	n/a	1


	205-3006-RLS			
	Adapter, 1/4" HB x 4			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	1/4"	4	1


	205-1008-RLS			
	Cap, Open top with adapter, 1/2" HB x 2 and Vent			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	1/2"	n/a	1


	205-3007-RLS			
	Adapter, 3/8" HB x 4			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	3/8"	4	1


	205-3001-RLS			
	Adapter, Closed			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	n/a	n/a	1

	205-5001-RLS			
	Adapter, 1/4" HB x 2 and Vent			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	1/4"	n/a	1

	205-3002-RLS			
	Adapter, 1/4" HB x 2			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	1/4"	2	1

	205-5002-RLS			
	Adapter, 1/2" HB x 2 and Vent			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	1/2"	n/a	1


	205-3003-RLS			
	Adapter, 3/8" HB x 2			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	3/8"	2	1


	205-3008-RLS			
	Adapter, Closed with Clear Window			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	n/a	n/a	1


HB = Hose Barb





VersaCap™ 120mm


	205-2003-RLS			
	Cap, Closed			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	120	n/a	n/a	2


	205-4005-RLS			
	Adapter, 3/8" HB x 3			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	120	3/8"	3	1


	205-2004-RLS			
	Cap, Open with closed adapter			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	120	n/a	n/a	2


	205-4006-RLS			
	Adapter, 1/4" HB x 4			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	120	1/4"	4	1


	205-4001-RLS			
	Adapter, Closed			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	120	n/a	n/a	1


	205-4007-RLS			
	Adapter, 3/8" HB x 4			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	120	3/8"	4	1

	205-4002-RLS			
	Adapter, 1/4" HB x 2			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	120	1/4"	2	1

	205-4008-RLS			
	Adapter, EZ-Top adapter			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	120	n/a	n/a	1

	205-4003-RLS			
	Adapter, 3/8" HB x 2			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	120	3/8"	2	1

	205-4009-RLS			
	Adapter, Closed w/Clear Window			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	120	n/a	n/a	1

	205-4004-RLS			
	Adapter, 1/4" HB x 3			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	120	1/4"	3	1

HB = Hose Barb



VersaBarb[™] and Spigot Technology



225-4001-RLS

VersaBarb[™], Spigot fitting, closed cap

Hose Barb Size	Thread Diameter	Qty/pk
n/a	1 1/8"	1



225-4005-RLS

VersaBarb[™], Spigot fitting, with 3/4" HB

Hose Barb Size	Thread Diameter	Qty/pk
3/4"	1 1/8"	1



225-4002-RLS

VersaBarb[™], Spigot fitting, with 1/4" HB

Hose Barb Size	Thread Diameter	Qty/pk
1/4"	1 1/8"	1



225-4006-RLS

VersaBarb[™], Spigot fitting, 3/4" Sanitary Connector

Hose Barb Size	Thread Diameter	Qty/pk
n/a	1 1/8"	1



225-4003-RLS

VersaBarb[™], Spigot fitting, with 3/8" HB

Hose Barb Size	Thread Diameter	Qty/pk
3/8"	1 1/8"	1



225-1001-RLS

Spigot, 2 pack

Hose Barb Size	Thread Diameter	Qty/pk
n/a	1 1/8"	2



225-4004-RLS

VersaBarb[™], Spigot fitting, with 1/2" HB

Hose Barb Size	Thread Diameter	Qty/pk
1/2"	1 1/8"	1

HB = Hose Barb



Chemical compatibility information is intended as a guide. We recommend that you test bottles for your specific application before use.

Roush Life Sciences, 8E Industrial Way #8, Salem, NH 03079 (603) 540-5778
www.RoushLifeSciences.com